

## pluriMate - 2 ml centrifuge tube



|                          | Order No.   | Size     | Prefilled |
|--------------------------|-------------|----------|-----------|
| <input type="checkbox"/> | 44-00002-15 | 500 pcs. |           |
| <input type="checkbox"/> | 44-00002-10 | 50 pcs.  |           |
| <input type="checkbox"/> | 44-00002-11 | 100 pcs. |           |
| <input type="checkbox"/> | 44-09102-15 | 500 pcs. |           |
| <input type="checkbox"/> | 44-09102-10 | 50 pcs.  |           |
| <input type="checkbox"/> | 44-09102-11 | 100 pcs. |           |
| <input type="checkbox"/> | 44-09302-15 | 500 pcs. |           |
| <input type="checkbox"/> | 44-09302-10 | 50 pcs.  |           |
| <input type="checkbox"/> | 44-09302-11 | 100 pcs. |           |
| <input type="checkbox"/> | 44-09202-15 | 500 pcs. |           |
| <input type="checkbox"/> | 44-09202-10 | 50 pcs.  |           |
| <input type="checkbox"/> | 44-09202-11 | 100 pcs. |           |

### Description


pluriMate was developed for optimal separation of leukocytes and peripheral blood mononuclear cells (PBMC) from whole blood and bone marrow. The key feature of pluriMate is the porous sponge incorporated at the bottom of the centrifuge tube. This barrier consists of high-grade polyurethane. It prevents you from time-consuming and laborious overlaying of the sample material. Anticoagulated blood or bone marrow can simply be poured directly from the blood sampling tube into the pluriMate tube. The porous barrier prevents mixture of the sample material with the separation medium. During centrifugation, leucocytes, lymphocytes and PBMCs are separated from unwanted erythrocytes and granulocytes, depending on the density gradient used (**Leuko Spin Medium**, **Lympho Spin Medium**, **Lympho 24+ Spin Medium** or **PLT Spin Medium**) on the basis of their density and enriched in an interphase above the separation medium. When separation is complete, the barrier prevents recontamination of the enriched cell fraction during harvest. Regardless of whether large or small sample volumes are to be processed. pluriMate is offered in 3 different tube sizes.

Use **only** with Whole Blood, Buffy Coat, Cord Blood, Bone Marrow.

## Application

Protocol for the use of pluriMate® centrifuge tubes

Use in combination with **density gradient media**.

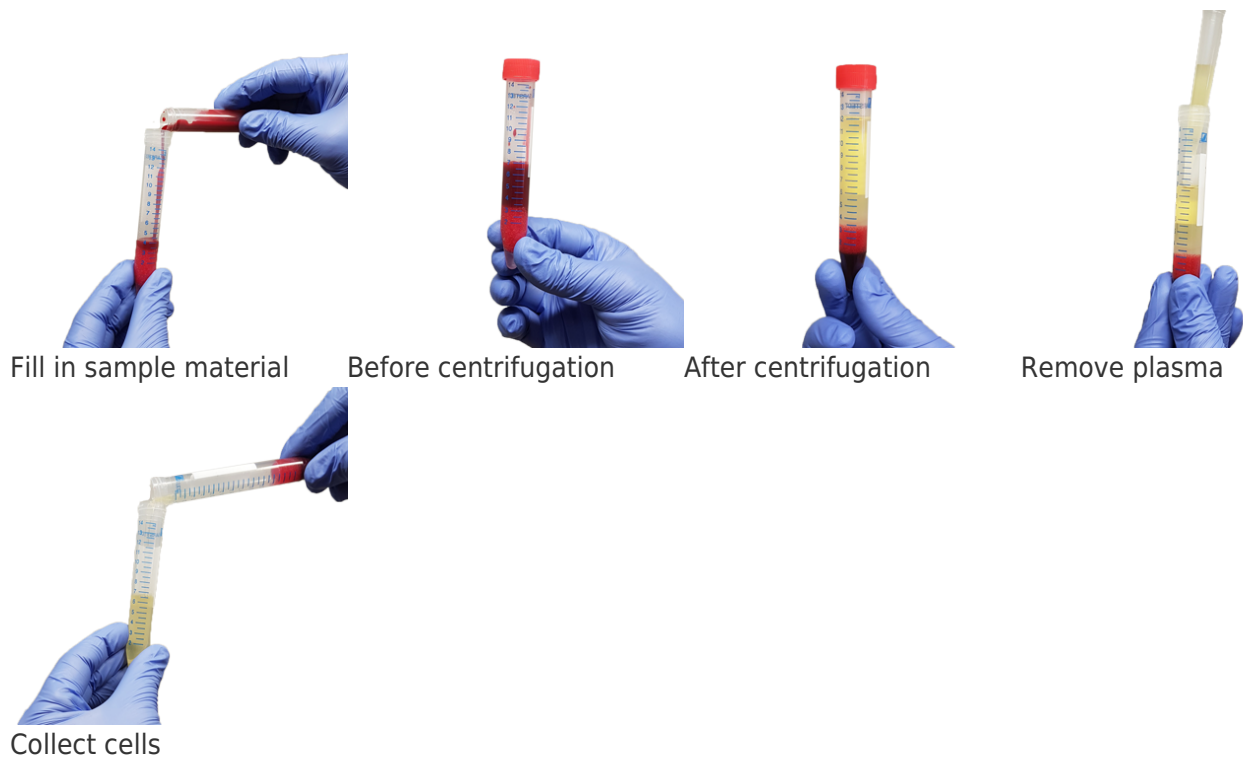

Find more information in the pluriMate® **protocol**.

## Additional Information

|                      |  |
|----------------------|--|
| Delivery Time (days) | 1-2  |
| Storage Condition    | Room Temperature   |
| Regulatory Statement | For research use only. Not for use in diagnostic procedures, |

## Warning and Limitations

This product is for research and development only, not for diagnostic or theurapeutic use.